



biomat 7000_{mx}



Read all precautions and instructions in this manual before using this Amethyst Biomat_{mx}. Keep this manual for future reference.



GENERAL INSTRUCTIONS FOR THE TMBIOMAT_{MX}

- 1 Insert White Connector from Controller Biomat_{mx} (QUEEN SIZE : L ▶ L / R ▶ R)
- 2 Insert electrical plug from controller into a wall outlet electrical socket.
- 3 Turn on the Control Box Safety Switch by pressing the button located on the right side of the unit.
- 4 Press Power Button to Activate the Biomat_{mx}. If you followed the above steps and the power is still not on, repeat steps 1-4 again.
- 5 Temperature Up / Down
Sets the desired temperature of the Biomat_{mx}. You will see the light move up or down on the indicator as you press the button.
- 6 TIMER :
Selects the length of time you want to use the Biomat_{mx}. The light will change to a different time setting as you press the Time Button on the indicator.
- 7 EMI(Electromagnetic Interception) :
Blocks harmful Electromagnetic and Water Vein Waves that would otherwise emit from the Biomat_{mx}.
This light will come and stays on.
- 8 TEMPERATURE INDICATOR :
This light will come on and stays on as you change the temperature control.
- 9 ALARM :
This light will appear if there is any malfunction. Please contact the dealer.

THERAPEUTIC RECOMMENDATIONS

DO NOT SLEEP ALL NIGHT USING A HIGH TEMPERATURE

- Children, diabetics, invalids or paralyzed persons should be monitored at all times while mat is in use. Never use high heat temperatures on children, or those who are pregnant or lactating. Avoid any heat on infants.
- Use a high temperature of (F 140j-158j / C 65j- 70j) for only 20~60minutes per-day to relieve arthritis pain and stress or for a relaxing sauna/detoxification.
- Use a low setting of (F 95j-113j/ C 35j- 45j) for a good nights sleep, or use no heat(power on only) for cool refreshing negative ions.



THE
RICHWAY
GROUP

Richway International Inc.
Richway Japan Co.
Richway EU Co.
The Life Corp.



EN ISO 13485 :
2003 Quality Management System
U.S. Trademark : 702502



biomat 7000_{mx}

BIOMAT_{MX} AND CONTROLLER

USER'S MANUAL

- Amethyst Biomat_{MX} Professional
- Amethyst Biomat_{MX} King
- Amethyst Biomat_{MX} Queen
- Amethyst Biomat_{MX} Single
- Amethyst & Tourmaline Biomat_{MX} Mini
- Amethyst & Tourmaline Pillow

SAFETY SYSTEM



RICHWAY INTERNATIONAL INC. www.richwayusa.com
1314 SOUTH KING STREET SUITE 520 HONOLULU HI 96814 USA. TEL: (808)589-2800 FAX: (808)597-1651



RICHWAY INTERNATIONAL JAPAN www.richwaynihon.co.jp
福岡県大野城市仲畑2丁目6-41 リッチウェイインターナショナル 日本支社 TEL: (092)404-3185 FAX: (092)404-3186



PRO BIOMAT & SINGLE BIOMAT (CONTROLLER)

NOTE ► PLEASE LOOK AT THE DRAWING, AND FAMILIARIZE YOURSELF WITH THE PARTS.

™ AMETHYST BIOMAT_{MX} PROFESSIONAL / ™ AMETHYST BIOMAT_{MX} SINGLE


RICHWAY INTERNATIONAL'S BIOMAT_{MX} CONTROLLER(SINGLE) FEATURES A HIGHLY DEVELOPED I.C CHIP WHICH CUTS OFF ELECTROMAGNETIC AND WATER VEIN WAVES, THE CONTROLLER IS DESIGNED WITH A SOFT TOUCH CONTROL SYSTEM, WHICH HAS MANY CONVENIENT FUNCTIONS.

Automatic overheating safety feature 2012 :

When temperature is set above 120°F / 49°C, Biomat will automatically reset itself down to 113°F / 45°C after 4 hours of use. This feature protects the user from overheating and dehydration.

FIRST,
INSERT THIS PLUG
INTO THE BIOMAT_{MX}

SECOND,
INSERT THIS
PLUG INTO
ELECTRICAL
SOCKET
ON THE WALL

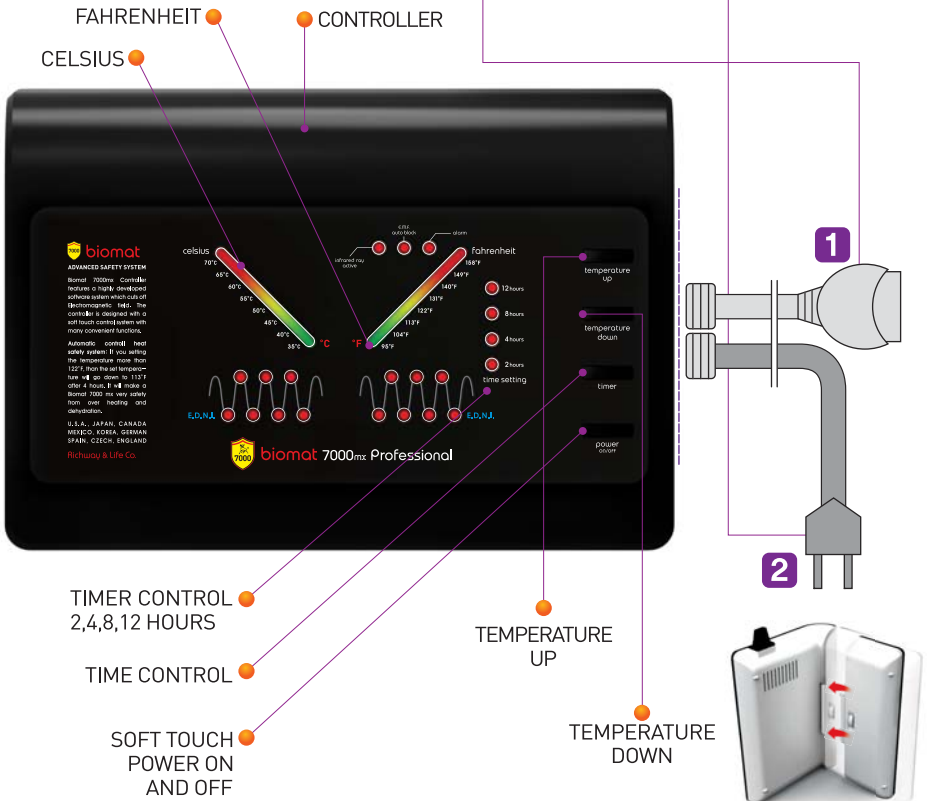



TM MINI BIOMAT_{MX} USER'S MANUAL


NOTE ▶ PLEASE LOOK AT THE DRAWING, AND FAMILIARIZE YOURSELF WITH THE PARTS.

AMETHYST & TOURMALINE MINI™ BIOMAT_{MX}

RICHWAY INTERNATIONAL'S™ MINI BIOMAT_{MX} CONTROLLER(SINGLE) FEATURES A HIGHLY DEVELOPED I.C CHIP WHICH CUTS OFF ELECTROMAGNETIC AND WATER VEIN WAVES.


FIRST,
INSERT THIS PLUG
INTO THE BIOMAT_{MX}

1

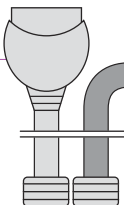

SECOND,
INSERT THIS PLUG
INTO ELECTRICAL SOCKET
ON THE WALL


2

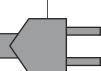

Automatic overheating safety feature 2012 :


When temperature is set above 120°F / 49°C, Biomat will automatically reset itself down to 113°F / 45°C after 4 hours of use. This feature protects the user from overheating and dehydration.

TEMPERATURE CONTROLLER
MINIMUM 104°F (40°C)
MAXIMUM 158°F (70°C)

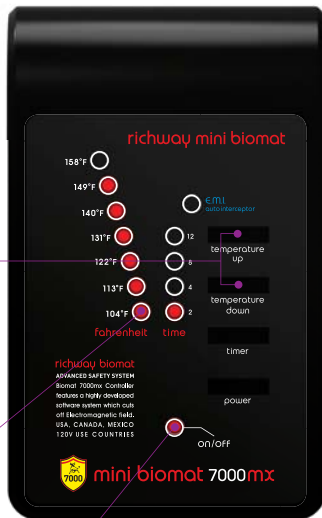

TEMPERATURE ON/OFF DISPLAY

POWER ON/OFF DISPLAY

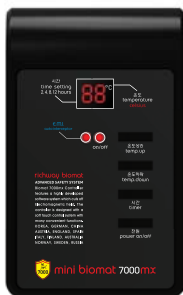

◀ 220 VOLTS MINI BIOMAT CONTROLLER (DIGITAL MONITOR TYPE)



Natural Infrared Rays and Negative Ions.

When you plug in electric outlet, if EMI light is not turn on, than switch the electric plug as left and right plug is changed.



QUEEN & KING BIOMAT (CONTROLLER)

NOTE ► PLEASE LOOK AT THE DRAWING, AND FAMILIARIZE YOURSELF WITH THE PARTS.



™ AMETHYST BIOMAT_{MX} QUEEN AND KING

RICHWAY INTERNATIONAL'S BIOMAT_{MX} CONTROLLER(QUEEN) FEATURES A HIGHLY DEVELOPED I.C CHIP WHICH CUTS OFF ELECTROMAGNETIC AND WATER VEIN WAVES, THE CONTROLLER IS DESIGNED WITH A SOFT TOUCH CONTROL SYSTEM, WHICH HAS MANY CONVENIENT FEATURES.

Automatic overheating safety feature 2012 :

When temperature is set above 120°F / 49°C, Biomat will automatically reset itself down to 113°F / 45°C after 4 hours of use. This feature protects the user from overheating and dehydration.

- FIRST,**
INSERT THIS PLUG INTO THE BIOMAT_{MX}
- SECOND,**
INSERT THIS PLUG INTO ELECTRICAL SOCKET ON THE WALL

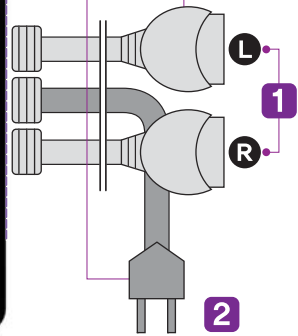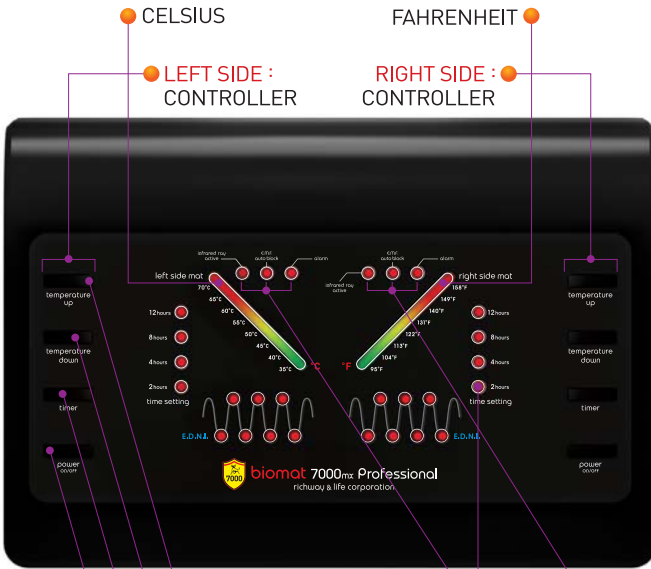

● CELSIUS FAHRENHEIT ●
● LEFT SIDE : CONTROLLER RIGHT SIDE : CONTROLLER ●

- LEFT SIDE : TEMPERATURE UP
- LEFT SIDE : TEMPERATURE DOWN
- LEFT SIDE : TIME CONTROL
- LEFT SIDE : SOFT TOUCH POWER ON AND OFF
- RIGHT SIDE : TEMPERATURE UP
- RIGHT SIDE : TEMPERATURE DOWN
- RIGHT SIDE : TIMER
- RIGHT SIDE : POWER ON/OFF
- RIGHT SIDE : ACTIVATE
- LEFT SIDE : ACTIVATE
- TIMER CONTROL 2,4,8,12 HOURS



PLEASE READ AND FOLLOW ALL OF THE SAFETY & OPERATING INSTRUCTIONS BEFORE USING YOUR BIOMAT_{MX}

▶ PLEASE READ FIRST

- 1 To reduce the risk of fire or electric shock, do not expose this product to rain or moisture.
- 2 If an alarm is sounded, or smoke is detected, immediately turn the power off and pull the power cord from outlet. Contact your sponsor and/or Corporate Office for assistance.
- 3 Do not position or suitable this product near heat sources or magnetic fields such as radiators, heat registers, stoves, TV sets, or magnetized objects.

▶ POWER SOURCE PROTECTION

Your Biomat_{MX} has a polarized plug (on blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug dose not fit fully in the outlet, reverse the plug. If it still dose not fit, do not continue attempting to plug in and contact a qualified electrician.

IF YOU NEED TO RETURN THE PRODUCT FOR ANY REASON, CONTACT CUSTOMER SERVICE FOR AN AUTHORIZATION CODE.


▶ REPLACEMENT PARTS

To order replacement parts, simply call our Customer Service Line at 1-808-589-2802, Monday through Friday, 9A.M. to 6P.M. Hawaii Time(excluding holidays). To help us assist you, please be prepared to give the following information when calling :

The MODEL NUMBER of the product (Single, Queen, BMP, or Mini)
The ORDER (RI) NUMBER. / Date of purchase / Applicant's Name

▶ SAFETY SYSTEM

- ▶ **TIMER** : To set the time want to use the Mattress. The light will change to a different time setting as you press the Timer Button on the indicator.
- ▶ **EMI (Electromagnetic Interception)** : Block the harmful Electromagnetic Waves.
- ▶ **TEMPERATURE INDICATOR** : When you change the temperature, this light will come on and stay on.
- ▶ **ALARM** : This light will appear if there is any malfunction. Please contact the dealer.

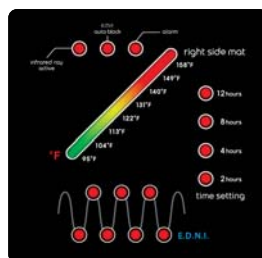

™BIOMAT_{MX} SAFETY INSTRUCTIONS- DO NOT DESTROY

(PLEASE KEEP THIS MANUAL IN A SAFE PLACE FOR FUTURE USE)
READ THE FOLLOWING IMPORTANT PRECAUTIONS BEFORE USING THE BIOMAT_{MX}

Read all instructions in this manual before using Biomat_{MX}. It is the responsibility of the owner to ensure that all users of the Biomat_{MX} are adequately informed of all precautions.

Place the Biomat_{MX} on a level surface, with a mat beneath it to protect the bed. Keep the Biomat_{MX} indoors, away from moisture and dust. Make sure plug is inserted firmly.

Replace any worn parts immediately. The Biomat_{MX} is intended for in home use only. Do not use Biomat_{MX} in commercial, rental, institutional setting. Keep children under the age 12, and pets away from the electric socket.



LIMITED WARRANTY

RichWay International, Inc. guarantees that for a period of one year from the date of purchase, this product shall be free of mechanical and electrical defects. Our obligation is limited to repair or replacement, at the manufacturer's option, of this product during the warranty period.

RichWay International, Inc. shall not be liable for any incidental or consequential damages as a direct result of negligence or misuse of this product, use on improper voltage/current, or by usage contrary to the operating instructions. Disassembly, repair or alteration, by any person other than those from authorized service center, is strictly prohibited. Mini and Pro mat's maximum charge is less \$180 in life time, Single, Queen and King's maximum repair charge is less \$240 in life time. Life time trade in service for all of Richway product.

▶▶ CAUTION

TO REDUCE THE RISK OF BURNS, ELECTRICAL SHOCK, FIRE AND ACCIDENT, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS.

- 1 READ ALL INSTRUCTIONS BEFORE USING THE MODEL RII-S1 / RII-Q2
- 2 It is the responsibility of the owner to ensure that all users of the Biomat_{MX} are adequately informed of all precautions.
- 3 Place the Biomat_{MX} on a level surface such as a bed or floor-mat to protect the product. Keep the Biomat_{MX} indoors and away from moisture and dust.
- 4 The Biomat_{MX} comes with a removable mattress pad. Do not use the mat without the padding in place. To ensure heat consistency, keep the Biomat_{MX} covered with a sheet, comforter or other natural fiber bedding at all times.
- 5 Do not use on areas of insensitive skin. Burns can occur regardless of control setting. Check skin frequently to avoid burning and blistering.
- 6 Do not use pins or other metallic means to fasten this pad in place.
- 7 Do not jump, stand, crush, or place heavy objects on pad. Avoid sharp folds and objects.
- 8 Never pull this pad by the electrical supply cord, control box or control box mat connectors. DO NOT USE ANY OF THE CORDS AS A HANDLE.
- 9 Do not use this pad with liniment, salve or ointment preparations that contain heat-producing ingredients. Skin burns could result.
- 10 Use this pad only with appropriate (110 or 220 volt AC circuit). Turn the pad off when not in use, or unplug when not in use over an extended period of time.
- 11 Do not tamper with this pad in any way. If for any reason this pad does not function satisfactorily, return it to the Corporate Office.
- 12 Do not immerse the pad or the controller in water. Clean with a damp cloth.
- 13 Do not use any detergents, solvents or flammable agents to clean the pad or the controller. DO NOT DRY CLEAN.
- 14 Unless under supervision, keep small children and pets away from the pad.
- 15 SAVE THESE INSTRUCTIONS.



AMETHYST PILLOW

NATURAL INFRARED RAYS

Pillow Size	480 × 300 × 120mm/ 19" × 12" × 4.3"
Net Weight	2.7kg / 6lb
Amethyst Weight	0.7kg
Tourmaline Weight	0.4kg

Elastic Polyurethane memory foam Pillow.

EXTERIOR DIMENSIONS AND WEIGHT

Electromagnetic Interception :
Auto

King Size	1800 x 2002 mm	72.65" x 77.81"	32kg	70.59 lb
Queen Size	1400 x 2002 mm	55.11" x 77.81"	23kg	50.70 lb
Single Size	900 x 2002 mm	35.43" x 77.81"	15kg	33.06 lb
BMP Size	720 x 1880 mm	28" x 74"	13kg	28.66 lb
Mini Size	508 x 812 mm	20" x 32"	3kg	6.61 lb

▶ RichWay International, Inc.™ BIOMAT^{MX} Dimensions

